

Non-invasive breath analysis for the detection of gastrointestinal disorders with the Gastrol^{yz}er[®] range



CE
0086



Helping to detect gastrointestinal disorders, one breath at a time.

www.bedfont.com



GastroCH₄ECK® Gastrolyzer® features

Accurate and real-time combined CH₄, H₂, and O₂ monitoring.





The role of Methane alongside Hydrogen Breath Testing (HMBT)

A percentage of patients do not produce hydrogen and therefore would not be diagnosed by the use of HBT alone. The role of GastroCH₄ECK® is to ensure that patients who are non-hydrogen producers, but produce methane or a combination of both, are not misdiagnosed.

“Medical literature shows that the level of methane producers varies by gender and population group, from approximately 33-41%^{5,6,7}”

How H₂ & CH₄ is produced

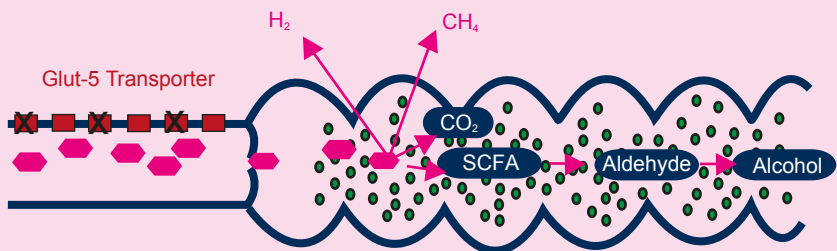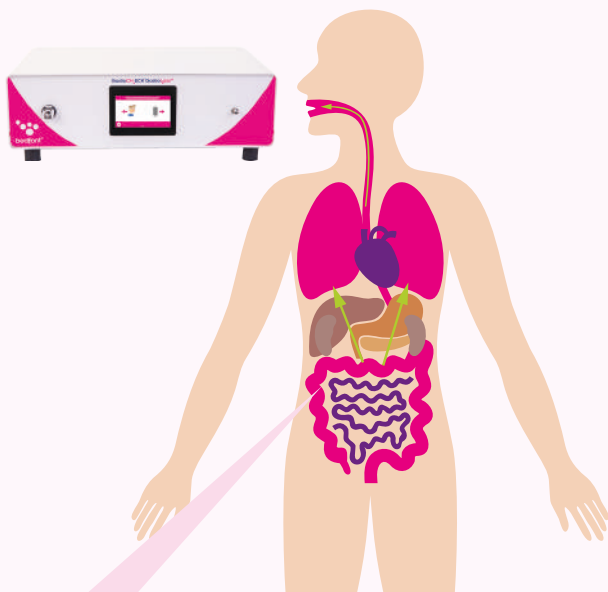

Consumables

SteriBreath™ mouthpieces: SteriBreath™ mouthpieces are both cost effective and compact. They are individually sealed for optimum infection control, whilst condensation in the tube confirms an accurate breath sample.



Order code: STERIBREATH-MP (250 per box)

D-piece™: The D-piece™ is used to attach a SteriBreath™ mouthpiece to the monitor. The D-piece™ incorporates a one-way valve and an infection control filter, which are proven to remove and trap >99% of airborne bacteria and >96% of viruses⁸. The D-piece™ should be changed every four weeks, an automatic reminder will appear on the screen every 28 days.



Order code: D-PIECE-3 (12 per box)

OneBreath™ mouthpieces: The OneBreath™ mouthpiece is a single-use bacterial filter mouthpiece and can be attached directly to the monitor to take a breath sample. The OneBreath™ incorporates a one-way valve and an infection control filter, which are proven to remove and trap >99.9% of airborne bacteria and <96% of viruses⁸.



Order code: ONEBREATH-MP (250 per box)

Face mask sampling system: This sampling system is single patent use and allows the patient to breathe normally through a face mask in order to produce a breath sample.



Order code: ISSA-V (1 sampling kit)

Infant face mask

Suitable for use with infants. **Order code:** EC60-IM-V

Medium face mask

Suitable for use with adolescent. **Order code:** EC60-MM-V

Adult face mask

Suitable for use with adults. **Order code:** EC60-AM-V



Monitor cleaning wipes: Free from alcohol to ensure continued performance of your monitor. Products containing alcohol cannot be used on any of the Gastrolyzer® range.



Order code: WIPE-V (pack of 50)

Calibration kit: The Gastro™ requires calibration quarterly using 100ppm H₂ gas, provided as a kit or replacement cylinder.

Order code: 012-14-12010K-V (kit)

012-14-12010-V (replacement cylinder)



GastroCH₄ECK® Gastrolyzer®

Hydrogen and methane monitoring to aid in the diagnosis of gastrointestinal disorders.

The GastroCH₄ECK® is a portable breath CH₄, H₂ and O₂ monitor, which allows health professionals to accurately detect a range of gastrointestinal disorders. An O₂ reading is taken to motivate patients to provide an end-tidal sample. If the sample is not adequate the GastroCH₄ECK® will automatically correct the reading saving the patient from embarrassment and eliminates the need to carry out another test.

Applications

The GastroCH₄ECK® can be used as an aid to diagnose the following disorders:

- Carbohydrate breakdown deficiency
- Carbohydrate malabsorption
- Lactose intolerance
- Bacterial overgrowth
- Determination of time passage through gut
- Irritable bowel syndrome (IBS)

Features/benefits

- Direct breath sampling line for instant results
- Breath bag sampling for simultaneous testing of larger groups
- Requires calibration just once a month, saving time and calibration gas costs
- Full-colour touchscreen
- Intuitive interface
- No warm-up time
- Maintenance reminders
- Service diagnostic management
- FREE GastroCHART™

“Medical literature shows that the level of methane producers varies by gender and population group from approximately 33-41%^{5,6,7}”

“If you are serious about providing a gold-standard breath testing service, the Bedfont® GastroCH₄ECK® system is essential, as without methane monitoring you will be missing very important clinical information and increase your percentage of false negative tests. It is also very important to map patients symptoms prior to and during the test to make sure the clinical translation of the physiological findings are maximised.”

Dr Anthony Hobson from The Functional Gut Clinic, London, UK.

GastroCH₄ECK® technical specification:

Concentration range	CH ₄	0-200ppm	
	H ₂	0-200ppm	
	O ₂	14-23%	
Power input		230V/100V, 50Hz-60Hz, 0.5-1.0 A	
Fuse		T 3.15 AH	
Calibration frequency		Every 4 weeks	
Start-up time		≤2 minutes	
Display		Full colour touchscreen	
Detection principle		Electrochemical sensor (O ₂ & H ₂) Laser (CH ₄)	
Accuracy	CH ₄	Resolution	1ppm
		Accuracy	±10% of reading
		Repeatability	<5% difference on consecutive readings
	H ₂	Resolution	1ppm
		Accuracy	±10% of reading
		Repatability	<5% difference on consecutive readings
	O ₂	Resolution	0.1%
		Accuracy	±10% of reading
		Repeatability	<5% difference on consecutive readings
Carbon monoxide cross-sensitivity (H ₂ only)		<4%	
Temperature range	Operating	15-35°C (59°-95°F)	
	Storage	0-40°C (32°-104°F)	
Pressure range	Operating	912-1114mbar (Atmospheric ±10%)	
	Storage	912-1114mbar (Atmospheric ±10%)	
Humidity range	Operating	30-75% RH (non-condensing)	
	Storage	15-90% RH (non-condensing)	
Sensor operating life	CH ₄	5 years	
	H ₂	2 years	
	O ₂	2 years	
Dimensions		Approx. 474 x 310 x 135 mm	
Weight		Approx. 8.5kg	
Monitor construction		Case: aluminium	
Classification		Class I ME equipment: (externally powered) Type BF applied part Method of sterilization (not suitable for sterilization) Not suitable for use in an oxygen rich environment Intended for continuous use	
Response time		≤45 seconds	

GastroCH₄ECK® consumables

GastroCH₄ECK® mouthpieces: specifically designed with the latest moisture removal and bacterial filtration to remove 99.9% of airborne bacteria from the patient's breath⁸. GastroCH₄ECK® mouthpieces are 'per patient use' so can be used for the duration of a patients testing procedure and disposed of. **DO NOT REUSE AFTER THE 1ST DAY OF TESTING.**



Order code: GASTROCHECK-MP-XL (box of 250)

GastroCH₄ECK® breath bags: specifically designed to capture as much end-tidal air as possible, with a one-way mouthpiece to prevent loss of sample. GastroCH₄ECK® breath bags are 'per patient use' so can be used for the duration of a patients testing procedure and disposed of. **DO NOT REUSE AFTER THE 1ST DAY OF TESTING.**



Order code: GASTROCHECK-BAG-XL (250 breath bags and bag mouthpieces)

Face mask sampling system: this sampling system is single patient use and allows the patient to breathe normally through a face mask in order to produce a breath sample. For use with breath bags.



Order code: ISSA-BB (1 sampling kit)

Infant face mask

Suitable for use with infants. **Order code: EC60-IM-V**

Medium face mask

Suitable for use with adolescent. **Order code: EC60-MM-V**

Adult face mask

Suitable for use with adults. **Order code: EC60-AM-V**



Moisture removal filters: for use with the breath bag sampling mode.

These will last for approximately 150 samples and need to be changed when the filter changes from orange to dark green.



Order code: GASTROCHECK-DESS-XL (pack of 5 moisture removal filters)

Calibration kit: the GastroCH₄ECK® requires calibration once per month using 100ppm H₂, 100ppm CH₄ & 20.9% air provided as a kit or replacement cylinder.

Order code: GASTROCHECK-CAL-KIT
GASTROCHECK-CAL -(Replacement cylinder)



GastroCHART™ patient database (included with every Gastro⁺™ & GastroCH₄ECK®)

The GastroCHART™ is designed specifically for use with the Gastro⁺™ and GastroCH₄ECK®. With the Gastro⁺™, the software is integrated and allows the health professionals to save up to 10 patients on the monitor for satellite breath tests. If the Gastro⁺™ or GastroCH₄ECK® are connected to a PC, readings can be downloaded to an unlimited patient database immediately. The readings can then be shown in tabular format or in a graph to show whether the patient has presented a positive/negative result. Results can be easily printed for the patient to retain and for record in their medical history.



References

1. Ledochowski, M. and Ledochowski, L. (2011) Hydrogen Breath Tests. 2nd Edition edn. Austria: Verlag Akademie für Ernährungsmedizin GmbH.
2. Eisenmann, A., Amann, A., Said, M., Datta, B. and Ledochowski, M. (2008) 'Implementation and interpretation of hydrogen breath tests', *Journal of Breath Research*, 2(4), p.046002.
3. Shelley, H., Brennan, M. and Heuschkel, R. (2009) 'Hydrogen breath testing in children: What is it and why is it performed?', *Gastrointestinal Nursing*, 7(5), pp. 18–27
4. Lee, W., Davidson, G., Moore, D. and Butler, R. (2000) 'Analysis of the breath hydrogen test for carbohydrate malabsorption: Validation of a pocket-sized breath test analyser', *Journal of Paediatrics and Child Health*, 36(4), pp. 340–342.
5. Roccarina, D., Lauritano, E. C., Gabrielli, M., Franceschi, F., Ojetti, V. and Gasbarrini, A. (2010) 'The role of methane in intestinal diseases', *The American Journal of Gastroenterology*, 105(6), pp. 1250–1256.
6. Pitt, P., de Bruijn, K. M., Beeching, M. F., Goldberg, E. and Blendis, L. M. (1980) 'Studies on breath methane: The effect of ethnic origins and lactulose', *Gut*, 21(11), pp. 951–954.
7. Di Stefano, M. and Corazza, G. R. (2009) 'Role of hydrogen and methane breath testing in gastrointestinal diseases', *Digestive and Liver Disease Supplements*, 3(2), pp. 40–43.
8. Public Health England. An Evaluation of Filtration Efficiencies Against Bacterial and Viral Aerosol Challenges Report No. 17/001. London: Public Health England; 2017

Join the **Gastrolyzer®** forum today

Did you know purchasing this product entitles you to free membership of the Gastrolyzer® email forum...

What is the Gastrolyzer® email forum?

It is an international, invitation only forum where professionals using the Bedfont® Gastrolyzer® range of monitors can communicate and share knowledge.

There is no cost or obligation to participate but membership is free when you purchase a Bedfont® Gastrolyzer® Gastro⁺™ or GastroCH₄ECK® - take advantage and join this elite group today!

For more information, please contact forums@bedfont.com.



How does it work?

Signing up is simple and only requires your full name and email address. We will then contact you for your Gastrolyzer® serial number.

Email forums@bedfont.com for more information.

Call on 01622 851122 or email ask@bedfont.com



Contact Bedfont® or one of our worldwide **Gastrolyzer®** distributors for a free demonstration.

www.bedfont.com
Tel: +44 (0)1622 851122
E-mail: ask@bedfont.com

A full list of our worldwide distributors can be found at
www.bedfont.com/distributors/

Our family, innovating health, for yours.

Bedfont® Scientific Ltd.

Station Road, Harrietsham, Maidstone, Kent,
ME17 1JA, England
Tel: +44 (0)1622 851122. Fax: +44 (0)1622 854860
Email: ask@bedfont.com
www.bedfont.com

© Bedfont® Scientific Limited 2019

Issue - 3 February 2019, Part No: MKT151_AX
Bedfont® Scientific Limited reserve the right to change or update this literature without prior notice.
Registered in: England and Wales. Registered No: 1289798

